


New!

For simplicity's sake

The Easy hygiene chair is designed to be simple, easy to use and easy to clean. Easy's user-friendly seat height mechanism allows adjustment to suit each individual user.

Easy is available in several models with different widths and with or without tilting. All models share the ability to use the newly developed comfort seat with the opening to the rear (accessory). The Easy (without tilting) can also be fitted with wheels (22") and hand rims to provide the user with the ability to move around on their own.

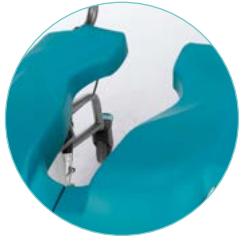

Integrated seat with the opening to the rear (accessory)

The seat has a soft, ergonomic design for the best comfort possible. The opening to the rear provides the carer with a natural, more hygienic working method.



Easy STD



Easy STD
Wheels and hand rims

Adjustable comfort back strap

Manufactured from a hygienic material – feels soft and comfortable. Easy to clean and never gets wet.

Tilt function

Easy Tilt can be operated either by gas strut or electric actuator (IP65).

Ergonomic seat

The Easy Tilt's seat is adapted for an upright posture that does not close the pelvis. It minimizes pressure and is designed to suit both men and women.

Plastic or metal

A guide is available for Easy as an accessory. Note that the guide is available in different versions to suit plastic or metal bedpans.

Adjustable seat height

An infinitely variable seat height from 48 all the way to 63 cm means the chair fits every patient and reduces staff workload to a minimum.

New height adjustment mechanism

Quick and easy to adjust seat height.

Durable material

Manufactured in corrosionproof materials.

Lockable castors

All four castors are lockable.

Swing-away, detachable footrests

Flip-up, swing-away, detachable footrests with heel restraints to prevent the patient's feet from sliding off in the tilt position are fitted as standard (Easy Tilt).



Technical specifications

CE marked according to MDD 93/42/EEC

Part No.	Product	Max user weight	Width/external dimensions	Internal dimensions/ armrests	Seat height	Seat depth	Tilt
56-270	Easy TILT Gas	150 kg	60 cm or 57 cm	49 cm or 44 cm	48-63 cm	45 cm	5 degrees forward to 25 degrees rearward
56-271	Easy TILT EI	150 kg	60 cm or 57 cm	49 cm or 44 cm	48-63 cm	45 cm	5 degrees forward to 25 degrees rearward
56-272	Easy STD	150 kg	60 cm or 57 cm	49 cm or 44 cm	48-63 cm	45 cm	N/A
56-273	Easy STD. Wheels and hand rims 22"	150 kg	75 cm	49 cm or 44 cm	49-64 cm	45 cm	N/A